



Safe hygiene and the fight against multi-resistant germs are becoming more and more important in care at home. Machine washable beds offer strong advantages here. They ensure **ideal cleaning and disinfection** with validated processes. At the same time, they save the operator time and money, since the time-consuming and arduous manual reprocessing by the staff of the medical supply company is no longer necessary. The bed can go back into circulation faster and unnecessary downtime is avoided.

The Dali wash was designed for cleaning in automatic washing systems. The bed can be reprocessed without disassembling the drive system because of its **IPX6 protection**. With its **CDP coating**, this bed is ideally protected inside and out against corrosion. CDP stands for "cathodic dip painting", an electrochemical painting process for maximum durability. Special drill holes allow residual water to drain off fully after the washing procedure.

Strong Advantages


- Safe working load 185 kg
- Full-length safety sides with Easy Click system
- Wireless handset
- Reverse-Trendelenburg included

Particularly high level of hygiene

Dali wash

The bed for optimal hygiene

Exterior dimensions/Height adjustment range:

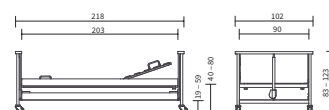

Product special features:

▲ The CDP coating protects the bed inside and out from corrosion.

▲ Drill holes in the metal surface allow any residual fluids to drain away.

▲ Plastic square plugs with drill holes allow residual moisture to escape.

▲ Sealed surfaces are ideal for safe and thorough cleaning.

▲ The 4-section mattress base has hygienic metal support bars.

