

ESSENTIALLY BETTER



Evario select

The ward bed for medical care

The Evario select is a model variant that enables good patient care at all times while focusing on the essential tasks in everyday hospital life. It is designed for safety, ease of use and the best possible support for the patient.

- Split safety side Protega with safe one-hand operation
- Integrated control panels unite the important functions for patients and nursing staff at a glance
- Safe working load of 225 kg
- Height adjustment range from 36 up to 84 cm depending on the type of castor
- Planar design for thorough manual reprocessing
- Attractive modern appearance with the 3 available head and footboard designs Intercontinental, Stelo and Rondo



- ▲ The control panels each have a simply designed user surface for the patient facing inwards.

Evario select

Safety and easy operation


The split safety side Protega allows either full-length protection or comfortable room for unhindered entry and exit. The Evario select is operated with panels integrated into the safety sides, with the options reduced to the most frequently used bed functions.

Advantages

- Protega safety side offers protection and promotes patient mobility at the same time
- Integrated panels with user surfaces located on the inner and outer sides enable simple and intuitive adjustment of the mattress base
- Nursing staff can trigger the CPR and Trendelenburg positions on the outer panels
- Optional control box for adjusting the bed from the foot end



- ▲ The series-standard control panels are located on both sides of the bed integrated into the Protega safety side.

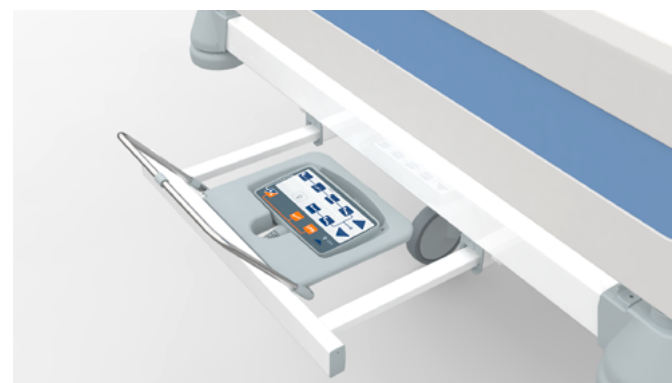

- ▲ Optionally the functions of the bed can also be operated with a pull-out control box at the foot end.



- ▲ The head-side elements of the safety sides move with the adjustment of the backrest. The patient is thus optimally protected and the control panels always remain within reach.

Evario select


Comfort and relief

Comfort for patients and relief for nursing staff are the strengths of the Evario select. Its height adjustment range from 36 up to 84 cm is ideal for a comfortable entry and exit as well as ergonomic work by the bed. The head-side elements of the Protega safety sides move along with the many health-promoting adjustment options of the mattress base and offer optimum protection. Tall patients benefit from the integrated bed extension.


Manoeuvring in confined spaces or a good straight line in long corridors are possible without any problems thanks to the robust 125 or 150 millimetre castors. An optional 5th castor further improves the handling characteristics. Colour-coded brake pedals facilitate the use of the castors' central locking mechanism.



- ▲ The integrated bed extension can be pulled out without tools.

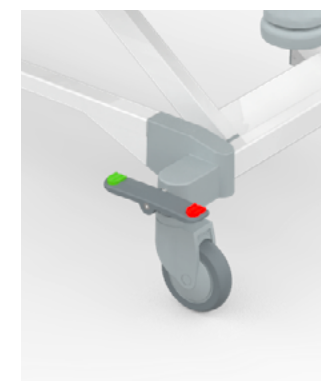

- ◀ The castors' central locking mechanism can be quickly activated or released with easily accessible brake pedals.

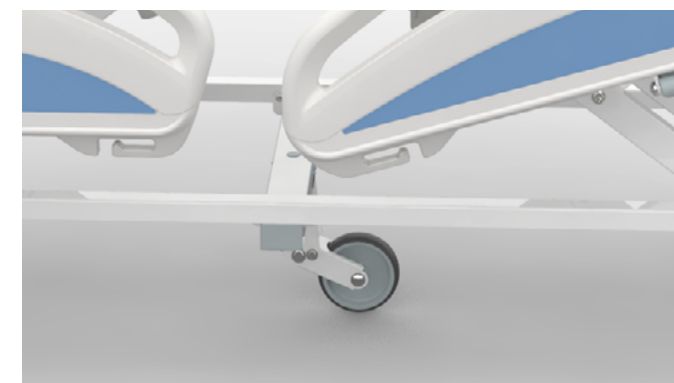

- ▲ The optional 5th castor further improves the good manoeuvring characteristics of the Evario select.

Equipment and technical data

Equipment / Options	
Safety sides	
Protega	✓
Head and footboard design	
Intercontinental	✓
Stelo	0
Rondo	0
Basic colour	
Crystal white	✓
Decontamination, degree of protection	
Manual reprocessing	✓
Mattress bases	
PP hygiene mattress base for belt system	✓
Wall deflection rollers	
4 plastic, all sides	✓
Castors	
Cover castor Ø 125 mm	✓
Integral castor Ø 150 mm	0
Brake pedal	
2 x plastic with coloured marking	✓
Directional castor position	
Head end, right	✓
Control units	
Integrated control panel ¹⁾	✓
Control box in the linen holder ¹⁾	0
Additional features	
5 th castor ¹⁾	0
Linen holder	0
Lower leg rest extension	0
Bed extension ¹⁾	0
Adapter sleeve for patient lifting pole	
Round pole	✓
Electronic features	
Foot pedal	0
Brake signal	0
Under bed light	0
Rechargeable battery	
Lead-gel battery	0

Technical data	
Dimensions and weights	
Overall dimensions	99.7 x 219 cm
Mattress base (mattress dimensions)	87 x 200 cm
Bed extension	28 cm
Mattress base division	83 – 15 – 33 – 53 cm
Backrest length compensation	10 cm
Safe working load	225 kg
Maximum weight of patient	ca. 185 - 215 kg
Ground clearance of chassis	ca. 22 - 25 cm
Empty weight	ca. 155 kg
Protective height of safety sides	41 cm
Adjustment ranges	
Tilting to Trendelenburg position	14°
Tilting to reverse-Trendelenburg position	14°
Mattress base height (Cover castor Ø 125 mm)	ca. 36 - 81.5 cm
Mattress base height (Integral castor Ø 150)	ca. 38.5 - 84 cm
Backrest angle	70°
Thigh rest angle	40°
Fowler position angle	16°



✓ = Standard 0 = Available as an option ¹⁾ = different type from the other Evario models