

ESSENTIALLY BETTER



Quado

The compact miracle

Always in the right place

The mobile bedside cabinet Quado

The new bedside cabinet Quado by Stiegmeyer is a real compact miracle. With its functions reduced to the most important ones it is particularly suitable for daily hospital routines. At the same time 2 open compartments on both sides and a drawer that can be opened to both sides offer enough storage space for the personal belongings of the patient. With its neutral colours Quado adapts well to room furnishings and every hospital bed. If the overbed table is not needed, it can be lowered easily and damped.



▲ With its neutral colours Quado adapts well to room furnishings and every hospital bed.




▲ The Quado's overbed table can be lowered easily and damped.



▲ The extendable bookrest under the overbed table top offers the patient high reading comfort, without having to clear the table (optional).

Transport

Bed with bedside cabinet

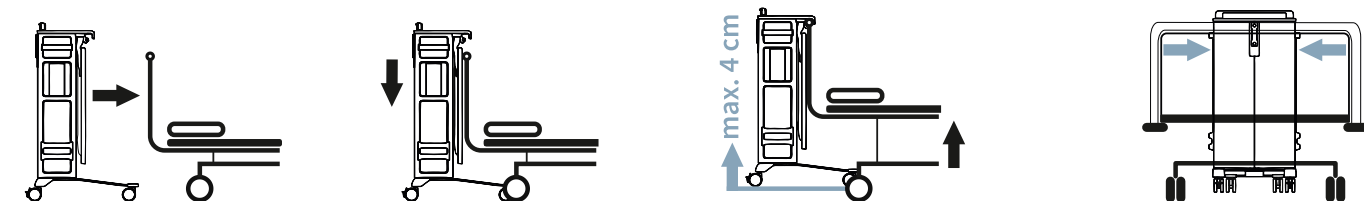

The easy transport of the Quado ensures relief during the daily routine. It can be attached to the head or footboard of almost all current Stiegmeyer hospital beds without effort, and thus can be transported together with the bed when a patient is relocated. The joint transport saves care staff from having to carry and, if applicable, document the personal items. During a simple transport the Quado can be rolled along with the bed. During a transport across thresholds or obstacles the cabinet should be lifted using the bed's

height adjustment. Once the Quado has arrived at its destination, it can be safely parked next to the bed thanks to its castors with brakes.



The Quado bedside cabinets ▶ can be pushed together for storage to save space.

Hygiene

Reliable reprocessing

The Quado is machine washable and can also be easily cleaned manually. This ensures optimum hygiene in the hospital.



▲ The chassis ensures high stability and straightforward cleaning.



▲ The hygiene-friendly design of the Quado lends itself to easy and thorough reprocessing.

Accessories

Useful additions



▲ With a towel hook on the rail of the bedside cabinet the towel is always at hand for the patient.

▲ With the longer hook toiletry bags can also be stowed close to the patient.

▲ An attachable USB extension cable allows charging of mobile devices in reaching distance of the patient.



► The drink holder offers secure hold for glasses and bottles.

◀ The optional additional compartment is suitable for magazines and further personal items.



▲► The optional sleeve for telephones and multimedia systems allows action in accordance with the patient's needs and the rooms' features.

Technical data / Equipment

Technical data		Castors	
Dimensions and weights		65 mm double castor (4 castors with brakes)	✓
External dimensions with overbed table (W x D x H)	92 x 44.6 x 89.6 –119.6 cm	Colours and materials	
External dimensions, overbed table retracted (W x D; incl. chassis)	56.4 x 44.6 cm	Body colour	Light grey NCS S1500-N
External dimensions body (W x D)	28.6 x 44.6 cm	Body material	PP, PA
Overbed tabletop height	85 –115 cm	Chassis colour	Chrystal White NCS S0500N
Overbed tabletop dimensions	40.7 x 57.4 cm	Chassis material	Steel
Total weight	ca. 25 kg	Miscellaneous	
Load capacity		Lock for top drawer	0
Overbed table (attached centrally)	max. 15 kg	Sleeves for telephone systems	0
Upper drawer pulled out (fixed centrally, tipping stability)	2,5 kg	Adaptor for extending the telephone sleeve	0
Top surface	15 kg	Suitable for automated cleaning and disinfection	✓
Equipment / Options		Bottle holder on top surface rail	0
Drawers		Towel hook on rail	0
Small top drawer	✓	Additional compartment on the side	0
Bookrest	0	Long hook on rail	0
		USB extension cable	0

✓ = standard feature
0 = available as an option

STIEGELMEYERGROUP



The Stieglmeyer Group is your international partner for hospitals and care facilities and, with the Burmeier brand, also for home care. We support our customers with practical products, solutions and services, which we constantly develop further.

- Founded in 1900 and still family-run today
- Employs 1,100 members of staff worldwide
- Production of all core products in our own factories
- Product development and manufacture in Germany
- Servicing, training and process optimisation at the customer's site
- International certification and export to more than 60 countries
- Sales companies in Belgium, Finland, France, the Netherlands, Poland and South Africa

Our brands





STIEGELMEYERGROUP

Hauptsitz / Headquarters

Deutschland / Germany
Stieglmeyer GmbH & Co. KG
Ackerstraße 42
32051 Herford
Phone +49 (0) 5221 185 - 0
Fax +49 (0) 5221 185 - 252
info@stieglmeyer.com
www.stieglmeyer.com

France
Stieglmeyer SAS
8 Allée des Chevreuils
69380 Lissieu
Phone +33 (0) 4 7238 - 5336
Fax +33 (0) 4 7238 - 5481
contact@stieglmeyer.fr
www.stieglmeyer.fr

South Africa
Stieglmeyer Africa (PTY) Ltd
The Estuaries · Century City
Cape Town, 7441
Phone +27 (0) 21 555 - 1415
Fax +27 (0) 21 555 - 0639
info.southafrica@stieglmeyer.com
www.stieglmeyer.co.za

Nederland / Netherlands
Stieglmeyer BV
Peppelenbos 7
6662 WB Elst
Phone +31 (0) 345 581 - 881
info@stieglmeyer.nl
www.stieglmeyer.nl

Suomi / Finland
Stieglmeyer Oy
Ahjokatu 4 A
15800 Lahti
Phone +358 (0) 3 7536 - 320
Fax +358 (0) 3 7536 - 330
stieglmeyer@stieglmeyer.fi
www.stieglmeyer.fi

Polska / Poland
Stieglmeyer Sp. z o.o.
86-212 Stolno
Grubno 63
Phone +48 56 6771 - 400
Fax +48 56 6771 - 410
info.stolno@stieglmeyer.pl
www.stieglmeyer.pl



Finden Sie Ihren Ansprechpartner auf unserer Website:
Find your local partner on our website:
www.stieglmeyer.com/partner

stieglmeyer.com