

Bedside cabinets



Clinic | Care | Homecare | Assist



Bedside cabinets
in a functional design





Contents



06 Hospital bedside cabinets



08 Vitano – variable and flexible



10 Optional equipment for Vitano



12 Somero – classic and stable



14 Optional equipment for Somero



16 Server Pleto – mobile and practical



18 Technical data and decors



19 Features and options



Hospital bedside cabinets

Even as a hospital patient, you still have all the little things with you that you need for everyday life. Your watch, your glasses, your smartphone – these are all indispensable personal belongings which need their own space. Stiegmeyer bedside cabinets are comfortably equipped so that you can reach your belongings easily from your bed. Drawers which only open at the front safeguard your property, while the closed rear panel protects your privacy.

These bedside cabinets are particularly stable, even if a heavy load is exerted on the table. They also support the work of nursing staff during many routine tasks.

Stiegmeyer bedside cabinets blend harmoniously with the design of hospital beds and complement other care furniture in the patient's room.



Vitano – variable and flexible

An ergonomic design is perfect for creating mobility. Take the Vitano bedside cabinet, for instance, with its ability to rotate. The body of the cabinet turns smoothly

and easily around a full 240°. Every patient can reach their bedside cabinet conveniently and safely. At the same time meals can be easily served and enjoyed.

Besides the obvious advantages, the Vitano bedside cabinet has many other useful details. These include features such as a towel hook, a rail to prevent objects

from falling off, and ergonomic handles. The well thought-out design of the chassis allows castors to be positioned with a high degree of stability.



Optional equipment for Vitano

Patients' needs and expectations are constantly expanding, and there is growing competition between hospitals. Attention is increasingly focused on the furnishings of rooms. The Vitano bedside cabinet, with its numerous versions, offers a variety of features and

provides comfortable storage space for the patient. It caters to individual wishes, and is available with features such as a fridge, a lockable drawer, a telephone holder or a fitting for modern multimedia monitors.



- Large picture: Fridge and lockable drawer
 Figure 1: Holders for phones and multimedia systems
 Figure 2: The overbed table is adjustable in height and can be pivoted. Overbed table and upper plate made of easy-care HPL in silver grey decor (optionally, same finish as front decor)
 Figure 3: Drink holder
 Figure 4: Double castors, optionally lockable
 Figure 5: Drawers with soft-close mechanism, with various options for storing personal belongings, additional drawer for slippers



Somero – classic and stable

The Somero bedside cabinet is designed to be suitable for use on either side of the bed. The overbed table can be moved from one side of the bedside cabinet to the other in a few easy steps due to its movable column. This offers complete flexibility in using the available space.

Somero also impresses with its sleek, clear-cut design. The column's additional castor makes the bedside cabinet easy to handle and enhances its stability. The Somero bedside cabinet is particularly strong and stable as a result.

- Large picture: The column and overbed table can be positioned on either side of the bedside cabinet
- Figure 1: Easy-care front through its sleek design
- Figure 2: Moving the column quickly and easily with one click of the foot pedal, without tools. The fifth castor makes it very stable.
- Figure 3: Towel hook as a useful detail



Optional equipment for Somero

Somero is a contemporary hospital bedside cabinet which can be extended in many ways. The optional

equipment provides practical solutions and additional features to meet a variety of needs and wishes.



1.



2.



3.



4.



5.



6.



7.

Large picture: Extra laptop drawer, optionally with a secure lock

Figure 1: Large fridge

Figure 2: The pivoting overbed table is height-adjustable, in the same finish as the front decor (standard)


Figure 3: Double castors, optionally lockable

Figure 4: Holders for phones and multimedia systems

Figure 5: Space for personal belongings

Figure 6: Slipper holder

Figure 7: Drink holder and towel hook

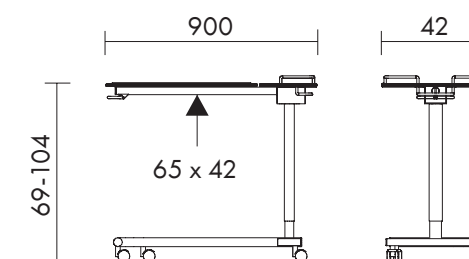

Server Pleto – mobile and practical

Due to its high mobility, the server Pleto can be universally used in daily nursing routines. Its main use is as a table. It can be positioned alongside beds, chairs and care chairs whenever necessary. Its advantageous C frame guarantees high strength and stability.

With its gas spring-assisted height adjustment, the server can be altered from heights of 69 cm to 104 cm. This allows it to be flexibly adapted to various uses, such as for eating or as an ergonomic bookrest

for reading. The table top can be horizontally pivoted. At the same time, a small section of the table remains horizontal, creating a practical surface on which objects can be placed. Patients can then, for example, read and put down a drink at the same time.

When the bed is raised, the server is automatically raised as well, since it can also be adjusted upwards without releasing the handle. Thus, collisions can effectively be prevented.


Technical data

Dimensions and weights	Vitano	Somero
External dimensions W x D x H	59.1 x 46.9 x 88.7 cm	55.5 x 44.4 x 87.3 cm
Height of overbed table	77 – 112 cm	76 – 110 cm
Dimensions of overbed table	HPL: 40 x 55.4 cm	HPL: 40 x 55.4 cm
Total weight	approx. 42 kg approx. 52 kg with fridge	approx. 38 kg approx. 46 kg with fridge
Load capacity		
Overbed table (fixed centrally)	max. 20 kg	max. 15 kg
Top drawer pulled out (fixed centrally, tipping stability)	max. 20 kg	max. 20 kg
Upper plate	max. 50 kg	max. 50 kg

Decor overview

To adapt the bedside cabinet to your own furnishing concept, we offer two body colours and a selection of different decors from the Resopal collection.

Please get in touch with us to let us help you plan your project.

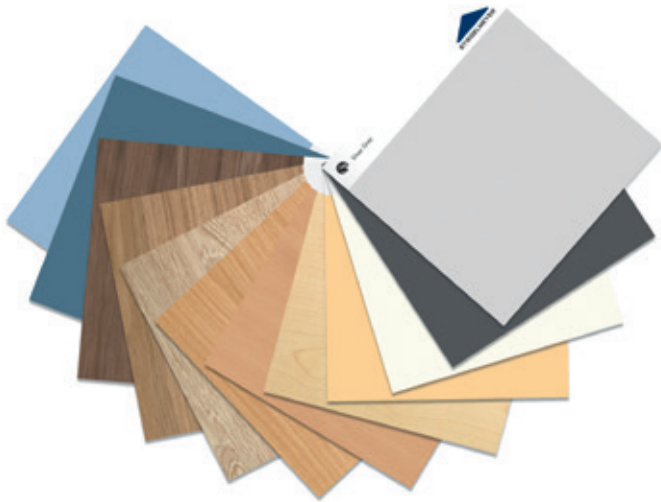

Features and options

	Model	Vitano	Somero
Drawers	Bottom drawer	●	●
	Small top drawer	●	●
	Large top drawer	–	○
	Large top drawer with integrated laptop drawer pullout	–	○
	Fridge	○	○
Colours and materials	Table and upper plate in HPL, silver grey decor	●	–
	Table and upper plate in HPL, with same finish as front decor	○	●
	Fronts, decorative laminate panels from Resopal collection	○	●
	Fronts, powder-coated aluminium	●	–
	Body colour argentum	○	○
Castors	Body colour crystal white	○	○
	52 mm double castor	–	●
	65 mm double castor	●	○
	2 castors with brakes	○	○
	4 castors with brakes	○	○
Positioning towards patient	Body can be rotated	●	–
	Overbed table column, movable	–	●
Miscellaneous	Lock for top drawer	○	○
	Slot-in holders for telephone systems	○	○
	Suitability for automatic cleaning and disinfection	●	○
	Bottle holder on upper plate rail	○	○
	Towel hook on rail	●	○
	Compartment for bedpan or slippers	○	●
	Additional side edges on upper plate	–	○
	Slipper holder	–	○

● = standard
○ = optional
– = not available

Ihr Bett – unsere Leidenschaft
Your bed – our passion

Finland / Suomi

Stieglmeyer Oy
Ahjokatu 4 B
15800 Lahti
Phone +358 (0) 3 7536 - 320
Fax +358 (0) 3 7536 - 330
E-Mail stieglmeyer@stieglmeyer.fi
finland.stieglmeyer.com

Netherlands / Nederland

Stieglmeyer BV
Peppelenbos 7
6662 WB Elst
Phone +31 (0) 345 581 - 881
E-Mail info@stieglmeyer.nl
netherlands.stieglmeyer.com

Poland / Polska

Stieglmeyer Sp. z o.o.
86-212 Stolno, Grubno 63
Phone +48 56 6771 - 400
Fax +48 56 6771 - 410
E-Mail info.stolno@stieglmeyer.pl
poland.stieglmeyer.com

France

Stieglmeyer SAS
451 Chemin de Champivost
69 760 Limonest
Phone +33 (0) 4 7238 - 5336
Fax +33 (0) 4 7238 - 5481
E-Mail contact@stieglmeyer.fr
france.stieglmeyer.com

South Africa

Stieglmeyer Africa (PTY) Ltd
The Estuaries, Century City, Cape Town, 7441
Phone +27 (0) 21 555 - 1415
Fax +27 (0) 21 555 - 0639
E-Mail info.southafrica@stieglmeyer.com
south-africa.stieglmeyer.com

Stieglmeyer GmbH & Co. KG

Ackerstraße 42, 32051 Herford, Germany
Phone +49 (0) 5221 185 - 0
Fax +49 (0) 5221 185 - 252
E-Mail info@stieglmeyer.com
www.stieglmeyer.com



Finden Sie Ihren Händler auf unserer Webseite:
Find your local partner at our website:
stieglmeyer.com/partner

Technische Änderungen vorbehalten. Leichte Farbabweichungen der Originale sind drucktechnisch bedingt. Stand 12/2016. Version 4. Artikel-Nr. 600576.
Subject to technical changes without notice. The shown product colours may differ slightly from the actual product. Last updated 12/2016. Version 4. Item no. 600576.